



E-Chrom Tech Co., Ltd.

TEL: +886-2-2511-1887

FAX: +886-2-2537-1870

EMAIL: service@echromtech.com


ADD: 10F-12, NO.20, Min Chuan W. Road, Taipei
104, Taiwan

WEB: www.echromtech.com

B10- Water Bath



WB-5000

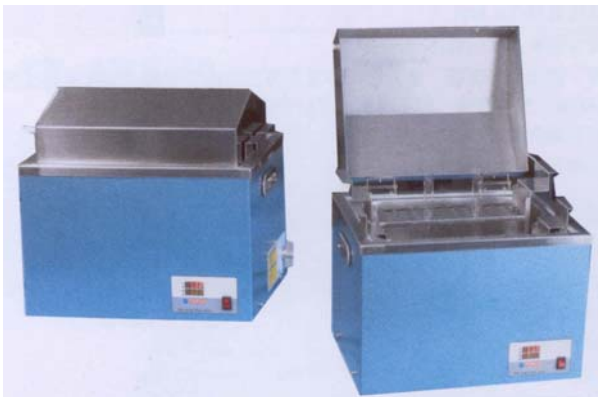

WB-10000



WB-20000

Descriptions

Water bath is an important thermostat equipment in chemical experiments. Because its high precision controllable heating circuit, the precision can reach to $\pm 0.2^{\circ}\text{C}$. Its LED display, on the other hand, provides clear and precise display of actual temperature and set temperature. In addition, this product uses durable material to ensure high productivity and long-term use.

Features

- A. Precise digital thermal control to ensure stability and precision.
- B. Two 4-digits LED display to show set temperature and actual temperature.
- C. Stainless steel tank and cover guarantee durability and cleanness.
- D. Automatic temperature protection alarm guarantees safe and reliable operation.
- E. Magnetic stirring under the stainless steel tank ensures even temperature.



E-Chrom Tech Co., Ltd.

TEL: +886-2-2511-1887

FAX: +886-2-2537-1870

EMAIL: service@echromtech.com

ADD: 10F-12, NO.20, Min Chuan W. Road, Taipei
104, Taiwan

WEB: www.echromtech.com

Specifications

Model	WB-5000	WB-10000	WB-20000
Capacity (L)	6L	10L	20L
Temperature Range	Room temperature + 5°C ~ 100°C		
Temperature Accuracy	±0.2°C	±0.2°C	±0.2°C
Power Consumption (W)	350W	550W	1100W
Internal Dimensions (W × L × D) (mm)	300 × 150 × 150	300 × 250 × 150	490 × 295 × 150
External Dimensions (W × L × D) (mm)	320 × 175 × 260	320 × 280 × 260	500 × 300 × 300